

**to remotely control the image analysis computer:
download team viewer www.teamviewer.us**

image analysis
computer

973 646 032

password
(this changes sometimes)

1yhc55


or just sit at the computer.

reserve the computer here:

<http://www.brownbearsw.com/cal/AnalysisComputer>

need your files organized like this:

- condition 1
 - date 1
 - cell 1
- do this semi-automatically by:
Saving movies of each channel separately in a folder
(one folder per date, no additional tif files in the folder)
Then open “makeCMEanalysisFolders_twocolor.m”
in Desktop/Track dynamics/ImageAnalysisPipeline/ and
click run.
- GFP/GFPmovie.tif
 - RFP/RFPmovie.tif

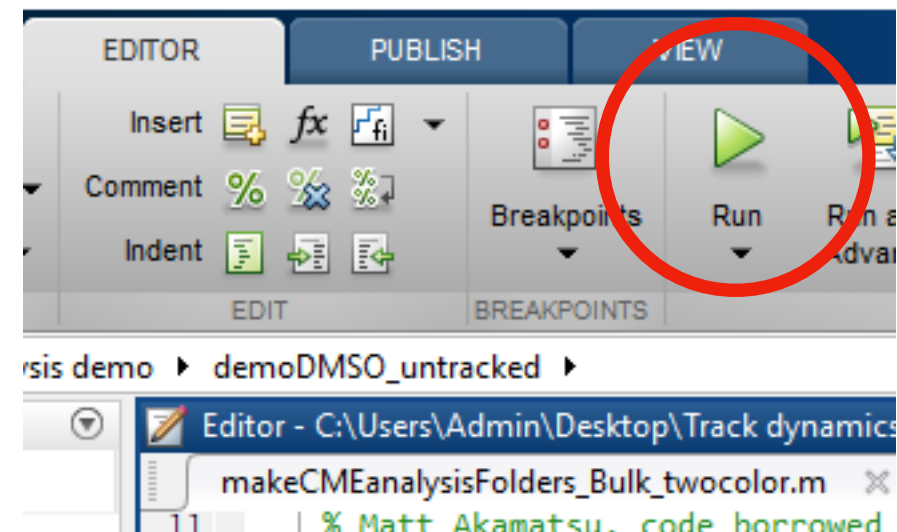


tracking

- open “runCMEanalysis.m”
- in Desktop/Track dynamics/ImageAnalysisPipeline/ .
click run.
- select the “condition” folder. (the “top” folder above
“date”). follow the prompts.
- Wait for it to complete. if you get an error you’ve
probably organized your files incorrectly.

plotting

- open up “plot_stats.m”
 - in Desktop/Track dynamics/ImageAnalysisPipeline/ . click Run.
 - click “pick analysis folder” to select your condition folder.
 - check “plot all cells”
 - type in time interval and PRESS ENTER
 - check the plot(s) you want.
 - all the “alignment” boxes plot the data with a different alignment, e.g. time 0 = when the reference track disappears.
 - “half dropoff” aligns to when one channel falls to 50% of its max intensity. I like using this alignment for dynamin.
 - click done.
 - if there’s an error your folders may not be organized properly.
 - you can save each graph with “file/save” or using the command `“saveCurGraphs(‘myExperiment’, [2 4 5])` and replace myExperiment with the name of your experiment, and the numbers with the figure numbers to save.

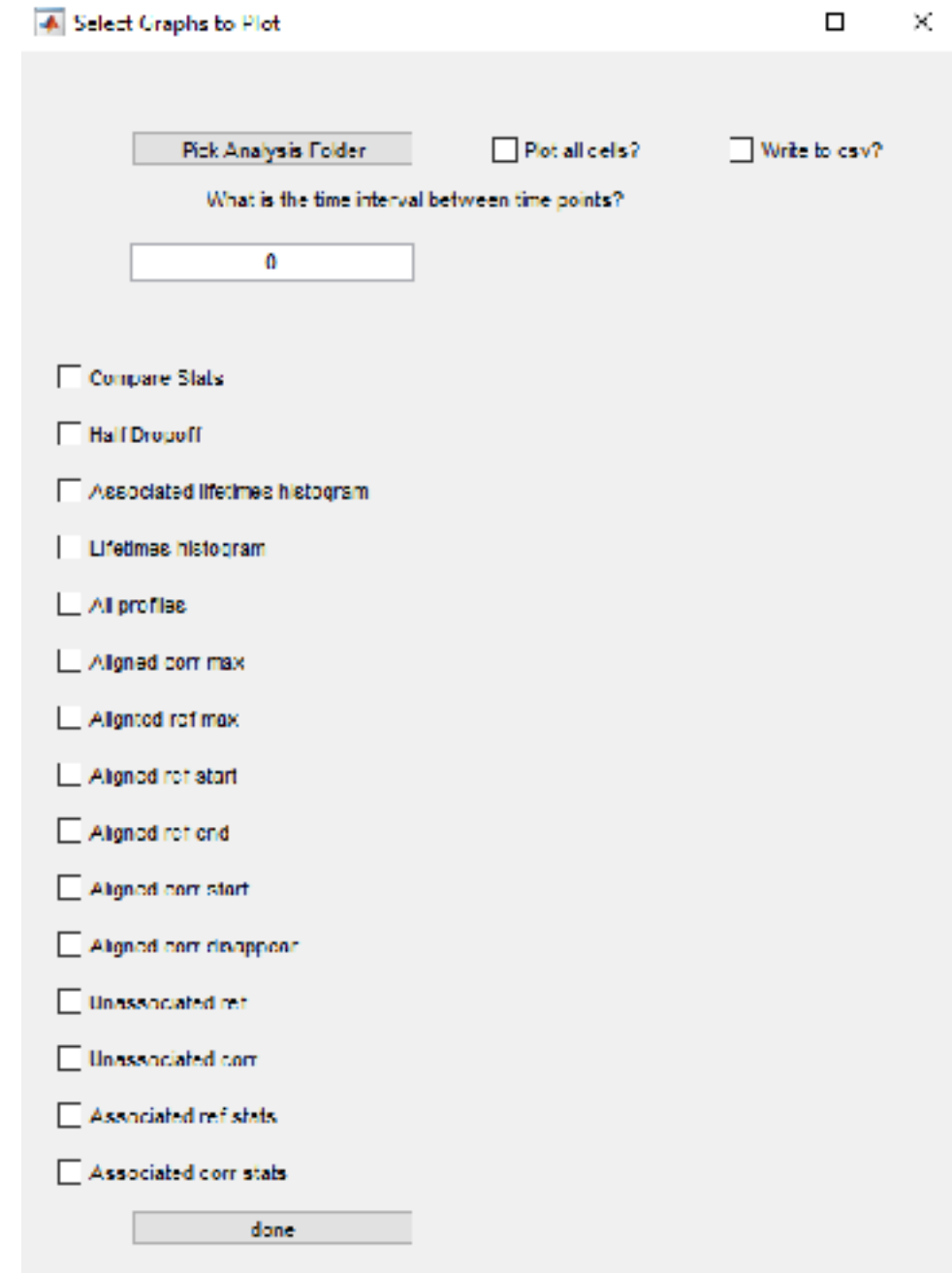

alternatively you can use the old version of plot_stats, in Desktop/Track dynamics/ImageAnalysisPipeline/AssociationCleaningPlotting

the csv writing and compare stats don't work yet

or: download the software to your own computer and use your own version of matlab

- download here: <https://github.com/DrubinBarnes/ImageAnalysisPipeline/tree/file-linking>

or clone the repository with this link:

<https://github.com/DrubinBarnes/ImageAnalysisPipeline.git>

password is yactin1

(this is unsupported)